



Oasis Series

# Y8y

## MULTIFUNCTIONAL ELECTRIC BED

Registration code: SKD-C




**SAIKANG** / Medical

Healthier  
your life

# Y8y

## FEATURES AND HIGHLIGHTS

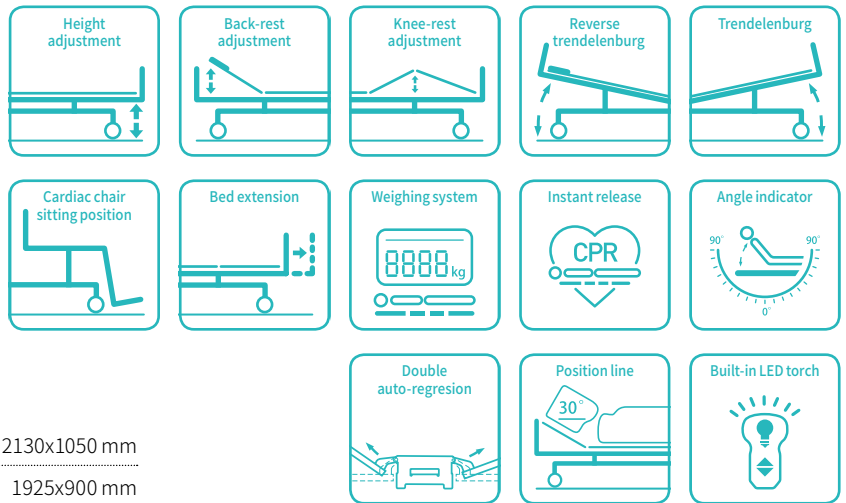

### Technical parameters

External size (LxW):	2130x1050 mm
Mattress platform:	1925x900 mm
Safe working load:	250 kg
Height adjustment:	350-800 mm
Back-rest adjustment:	0-75° (±10°)
Knee-rest adjustment:	0-35° (±10°)
Reverse Trendelenburg:	0-12°
Trendelenburg:	0-12°
Angle indicator:	0-90°
Bed Extension:	200mm
Auto-regression:	100+60mm
Position Line:	30°
Voltage and frequency:	100-240 V · 50/60 Hz
Plug as required (hospital grade)	

### Advantages

- ★ Double Auto-regression
- ★ Super-low Bed Position 350mm
- ★ PP Bed Platform, Removable & Washable
- ★ Inserted-type PP Head Foot Board, One Step Release With Safe Lock, Patent Design
- ★ PP Safety Guardrail With Two Steps Unlocking System, Patent Design
- ★ Back Rest 30degree Position, Reduce The Chance Of Lung Infection
- ★ Built-in LED Torch
- ★ Applicable For C-arm Examination

### Technical configuration


<input checked="" type="checkbox"/> Linear Motor	4pcs
<input checked="" type="checkbox"/> 5" Double Side Castors, Central Locking System	4pcs
<input checked="" type="checkbox"/> ABS Oasis Side Rail Extra Height	1set
<input checked="" type="checkbox"/> ABS Oasis Bed Ends with Safe Lock	1set
<input checked="" type="checkbox"/> PP Bed Platform, Removable & Washable	9pcs
<input checked="" type="checkbox"/> Plastic Socket Inside PP Guardrail for Phone or Glasses	2pcs
<input checked="" type="checkbox"/> Patent Angle Indicator For Back-rest And Trendelenburg Position	4pcs
<input checked="" type="checkbox"/> Ø100mm Wall Bumper Wheel	4pcs
<input checked="" type="checkbox"/> IV Pole	1set
<input checked="" type="checkbox"/> IV Pole Prevision	4pcs
<input checked="" type="checkbox"/> Drainage Hook	2pcs
<input checked="" type="checkbox"/> Central Locking Pedal	2pcs
<input checked="" type="checkbox"/> ABS Mattress Stopper	2pcs
<input checked="" type="checkbox"/> Hand Controller	1pc
<input checked="" type="checkbox"/> Manual CPR	

### Optional

<input type="checkbox"/> Electric CPR	
<input type="checkbox"/> Weighing System	
<input type="checkbox"/> Bed Extension 20cm	
<input type="checkbox"/> Back-up Battery	
<input type="checkbox"/> Nurse Control Panel	
<input type="checkbox"/> Emergency Stop Button	
<input type="checkbox"/> Cardiac Chair Position	
<input type="checkbox"/> Directional Wheels Function	
<input type="checkbox"/> 5th and 6th Extra Side Rail For Full Protection	
<input type="checkbox"/> Monkey Pole	
<input type="checkbox"/> Balkan frame	
<input type="checkbox"/> Linen Holder	

**Great Multi-Care**

Respects easy patient mobilisation: One-button for patient mobilisation position -back plate raised, double hole design on back part guardrail, which supports patient self standing.

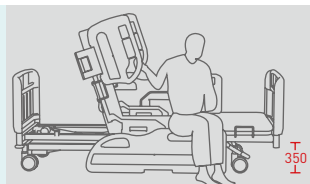


**Safety Guardrail**

The designed guardrail fully considered to provide enough space and hand grasp for patient to get out of bed, and provide extremely high safety. Two steps unlocking system, to prevent patients from accidentally lowering the guardrail (patented design).



**Double Regression Function**

Auto-Regression 10+6 cm  
When repositioning the backrest and thigh rest, the Y8Y8C makes the back plate retreat 10cm and the thigh plate retreats 6cm so to minimise the pressure and shear forces in the patient's pelvic region. This results in a significant reduction in the development of pressure ulcers.

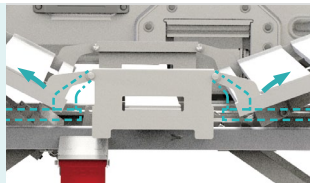

**Electric Motor**

Linear motor, with UL/ROHS/EN standard. Noise under 50 dB, with lifetime 20,000 times movement.



**Safety**

IEC60601 Safety Standard.  
One-button emergency stop switch can quickly stop all going movement.



**CPR**

CPR emergency reset handle to immediately flatten the bed back part; located on both sides of the bed. Electric CPR button is setted inside control panel, one button to make the back plate flat.



**Wheels**

Castors with TPR tire, no worn-out after running 30KM, united forming without bolts  
Pass Dynamic test: Bearing 120kg run 30KM and pass obstacles 500time.




**Brake**

Three position braking system, transporting can maintain a straight line.



**Joint-Free Technology**

Removable plastic covers formed from 9 removable blow moulded lamella to minimise the slots where unwanted infections may hide.

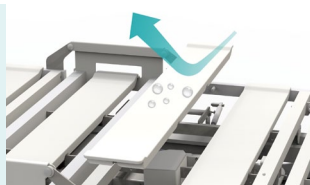

**Angle Indicator**

Patent angle indicator, clearly show back-rest and tredelenbug position.




**Head & Foot Board**

Patented design Detachable head & foot board, using the eccentric locking structure, one step quick removable.



**Bed Extension**

Extended bed platform up to 20CM, comfortable for the tallest of patients (option)

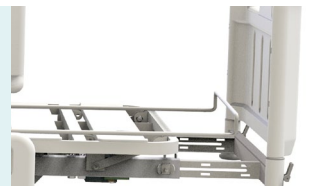

**Mattress**

Mattress Holder.



**Drainage Hooks**

Horizontal corner bumpers at the outermost periphery of the bed to protect the bed and wall.



**IV Pole**

Full stainless steel IV pole, bearing 15 kg.



**Drainage Hooks**

Mobile plastic drainage hooks.






## CONSTANT INNOVATION



Design



Research



Manufacturing



Maintenance



Service



## WORLDWIDE BRANCHES

### USA

**Saikang Medical Technology, Inc.**  
17800 Castleton St Ste 180,  
City of Industry, 91748, California.  
Tel.: +1 949648 2776

### United Kingdom

**UK Healthward International Ltd.**  
Dept. 111, 196 High Road Wood  
Green, N22 8HH, London.  
Tel.: +44 2032399738

### Russia

**Saikang Medical Russia LLC**  
Truda 174,  
454080 Chelyabinsk.  
Tel.: +7 (351) 214 5566 / 8811 / 1516

### East Africa

**Saikang Medical Solutions Ltd.**  
Viraj Complex Unit 3, Mombasa  
Road, P.O. Box 14262-00800, Nairobi.  
Kenya.  
Tel.: +254 733 704 690

### Middle East

**Saikang Medical Supplies LLC**  
X22-S03 England Cluster,  
International City, Dubai,  
United Arab Emirates.  
Tel.: +20 128 8962497

# SAIKANG / Medical

Jiangsu Saikang Medical Equipment Co., Ltd.  
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China  
 📞 +86 18021231901 · 🕒 24/7  
 ✉️ sales@chinasaikang.com  
 🌐 www.saikangmedical.com